

About Samsung India Electronics Pvt. Ltd.  
Samsung India Electronics Pvt. Ltd. inspires the world and shapes the future with transformative ideas and technologies, redefining the worlds of TVs, smartphones, wearable devices, tablets, cameras, digital appliances, medical equipment, network systems and semiconductors. We are also one of the lead in the Internet of Things space through, among others, our Digital Health and Smart Home initiatives. We employ 307,000 people across 84 countries. To discover more, please visit our official website at [www.samsung.com](http://www.samsung.com).

Copyright© 2025 Samsung India. All rights reserved.

Images of the foetus are digitally created for representational purposes. Samsung does not encourage prenatal sex determination, and the same is illegal in India. Image simulated for representational purposes only.

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer. Please contact your Samsung representative for the most up-to-date product information.

The information provided in this brochure is intended for general informational purposes only. Samsung disclaims all warranties, express or implied, related to merchantability, fitness, output of device, for any particular purpose or usage.

\*This product, features, options, and transducers may not be commercially available in some countries.

\*Sales and Shipments are effective only after the approval by the regulatory affairs.  
Please contact your local sales representative for further details.

\*This product is a medical device, please read the user manual carefully before use.

\*S-Vue Transducer™ is the name of Samsung's advanced transducer technology.

\*Optional feature which may require additional purchase.

# SAMSUNG

Relentless Innovation  
for your diagnostic confidence

## cV5

## For Cardiovascular Care



Contact us



WhatsApp



Samsung India Electronics Pvt. Ltd.

CIN: U37900DL1995PTC07138

Regd. Office: 20th-24th Floor, Two Horizon Center, Golf Course Road,

Sector-43, DLF Phase-5, Gurgaon-122002, Haryana (India)

Visit: [www.samsunghealthcare.com](http://www.samsunghealthcare.com) | E-mail: [time.marketing@samsung.com](mailto:time.marketing@samsung.com)

# Excellence in Cardiovascular Care, Tailored to Your Needs

Discover a new frontier in cardiovascular care with Samsung cV ultrasound systems. Our advanced tools are designed to empower clinicians with precise, high-resolution imaging, streamlining diagnostics and inspiring confidence through the latest advancements.

With a dedicated cardiovascular control panel, the system offers an intuitive, user-friendly experience, while a 68.58cm OLED monitor delivers enhanced clarity for critical insights.



## Key Benefits



### Outstanding image quality for confident diagnoses

Advanced imaging technology, characterised by its clarity and high resolution, enhances the accuracy of cardiovascular assessments.



### Comprehensive suite of advanced diagnostic tools

A versatile range of solutions, driven by automation and AI, helps elevate clinical confidence and improve efficiency.




### Tailored for seamless cardiovascular workflow

A system designed to meet the specific requirements of cardiovascular operations, ensuring optimal performance.




# Outstanding Image Quality for Confident Diagnoses


Advanced imaging technology, characterised by its clarity and high resolution, enhances the accuracy of cardiovascular assessments.

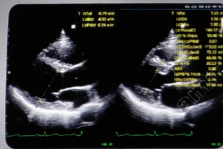

LV Study with  
HeartAssist™

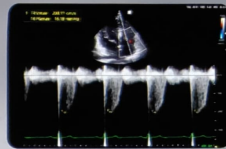

TR Average Velocity in  
CW Doppler with HeartAssist™

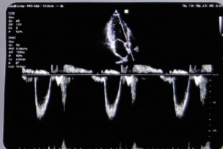

LVOT with  
PW Doppler



Carotid Artery with  
LumiFlow™




## Comprehensive Suite of Advanced Diagnostic Tools

A versatile range of solutions, driven by automation and AI, helps elevate clinical confidence and improve efficiency.




### Strain<sup>+</sup>


Strain<sup>+</sup> is a quantitative tool for measuring the global and segmental wall motion of the left ventricle (LV), left atrium (LA) and right ventricle (RV) to assess systolic and diastolic heart function.

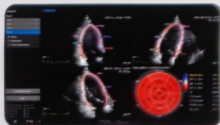

Strain+ LV

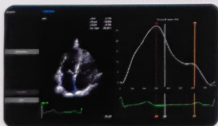


Strain+ LA

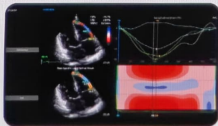

Strain+ RV


### HeartAssist<sup>™1</sup>

HeartAssist<sup>™</sup>, a feature that uses deep learning technology, automatically classifies ultrasound images into the required measurement views for heart diagnosis and provides measurement results.

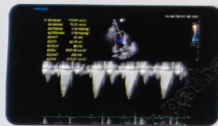

### AutoEF

AutoEF is a feature that conveniently measures and quantifies Ejection Fraction, LV volume and Global Longitudinal Strain (GLS). It calculates end-systolic and end-diastolic points of the left ventricle to facilitate efficient assessment of heart function.

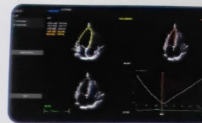


### StressEcho<sup>1</sup>

The StressEcho package features wall motion scoring and reporting. It includes exercise StressEcho, pharmacologic StressEcho, diastolic StressEcho and programmable StressEcho.




### AutoIMT<sup>+</sup>

AutoIMT<sup>+</sup> is a screening tool that analyses a patient's potential risk for cardiovascular disease. It enables easy measurement of intima-media thickness for both the anterior and posterior walls of the common carotid at the click of a button.

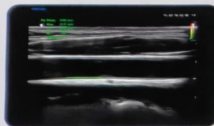

### ArterialAnalysis<sup>™</sup>


ArterialAnalysis<sup>™</sup> detects functional changes of vessels, providing measurement of stiffness, intima-media thickness, pulse wave velocity of the common carotid artery. Since functional changes often precede morphological changes, this technology enables the early detection of cardiovascular disease.

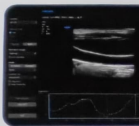




## Tailored for Seamless Cardiovascular Workflow

A system designed to meet the specific requirements of cardiovascular operations, ensuring optimal performance.

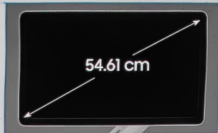

54.61 cm LCD Monitor\*

\*Monitor size/type may vary for cV8/cV7/cV6/cV5 ultrasound system.  
(68.58cm OLED/60.45 cm LCD/54.61 cm LCD)

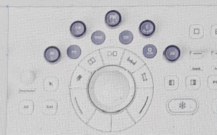

Easy-to-access CV modes with intuitively designed control panel

### Optimise image with just one click

QuickScan™ technology intuitively optimises grey scale and Doppler parameters. It allows users to adjust the ROI box location with a single press of a button.



### Choose transducer and preset combinations with a single click

Users can select the most common transducer and preset combinations with one touch. QuickPreset increases efficiency, making a full day of scanning simple and easy.



### Assign specific functions to the buttons near the trackball

Depending on the ultrasound inspection items, the functions assigned to the buttons surrounding the trackball can be utilised to minimise the hassle of menu selection.



## Comprehensive Selection of Transducers

### Phased Array Transducers



PA1-5A  
Abdomen, Cardiac,  
Paediatric, Thoracic,  
Vascular, TCD\*



PA2-9S  
Abdomen, Cardiac,  
Paediatric, Vascular

### Curved Array Transducers



CA1-7SD  
Abdomen, MSK, Vascular,  
OB, GYN, Thoracic,  
Paediatric, Urology



CA4-10M  
Abdomen, Vascular,  
Paediatric

### TEE Transducers



HMPT3-7  
Cardiac



TA2-9  
Cardiac

### Linear Array Transducers



LA2-9S  
Abdomen, MSK,  
Small Parts, Vascular,  
Paediatric



LA3-14AD  
Abdomen, MSK,  
Small Parts, Vascular,  
Paediatric



LA3-22AI  
MSK, Intraoperative



L3-22  
MSK, Small Parts,  
Vascular, Paediatric

### CW Transducers



CW6.0  
Cardiac, Vascular,  
TCD



DP2B  
Cardiac, Vascular,  
TCD

\* Supported transducers may vary depending on the cV8, cV7, cV6, and cV5 ultrasound system.

\* For detailed information regarding transducer compatibility and support for this product, please contact your sales representative.