

Relentless Innovation
for your diagnostic confidence

SAMSUNG



V6

Inspiring your everyday



Product Inquiry

Embracing efficiency in your daily ultrasound scanning

Begin your journey towards efficient healthcare with the Samsung V6 ultrasound system. Our robust solution for general imaging offers both image clarity and advanced automated features. Additionally, Samsung's cutting-edge imaging engine, Crystal Architecture™ ensures a reliable ultrasound experience.

Experience simplicity with our easy-to-use system, specifically designed to alleviate your workload and enhance usability. Furthermore, our powerful system comes with battery capability, providing additional operational convenience. The Samsung V6 ultrasound system is a partner you can depend on to deliver exceptional efficiency to meet your daily ultrasound needs.






View webpage





Elevating confidence with superb imaging performance

The V6 delivers exceptional 2D and color image quality tailored for general imaging, driven by Samsung's core imaging engine, Crystal Architecture™. With its comprehensive imaging capabilities, the V6 is designed to seamlessly support your daily ultrasound scanning needs, enabling clear and accurate image acquisition. Experience confidence and accuracy in ultrasound scanning with the V6.

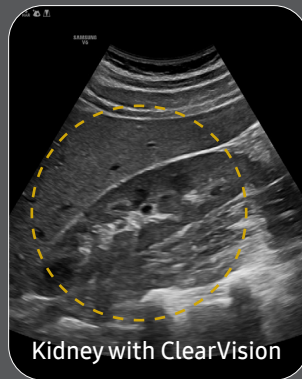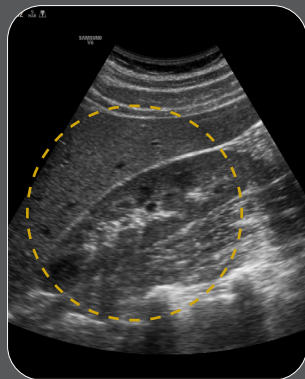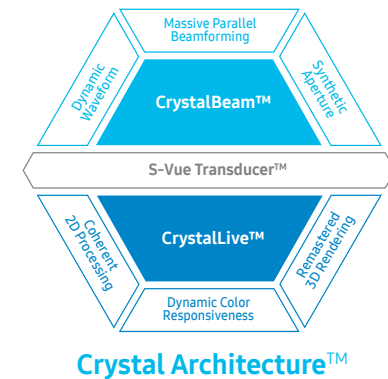



Kidney with ClearVision

Reduce noise to improve 2D image quality

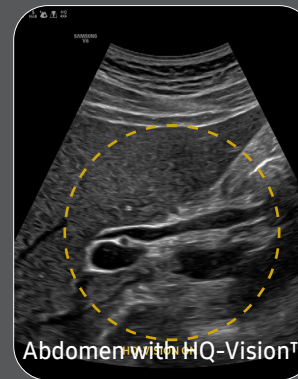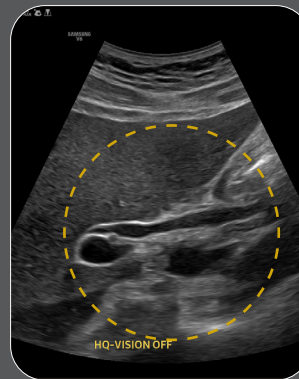


Abdomen with HQ-Vision™

Clean up blurry areas in the image

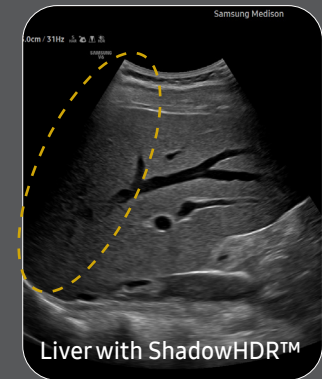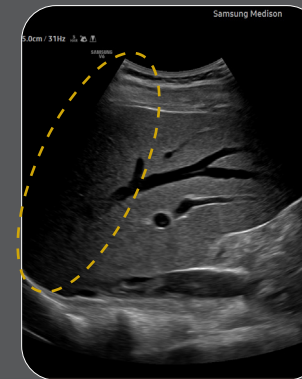


Liver with ShadowHDR™

Enhance hidden structures in shadowed regions

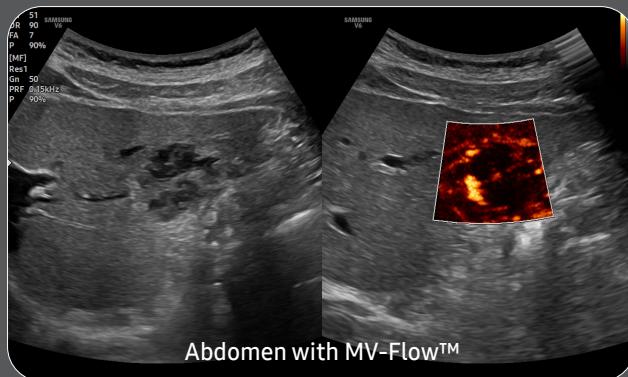


Abdomen with MV-Flow™

Visualize slow flow in microvascular vessels

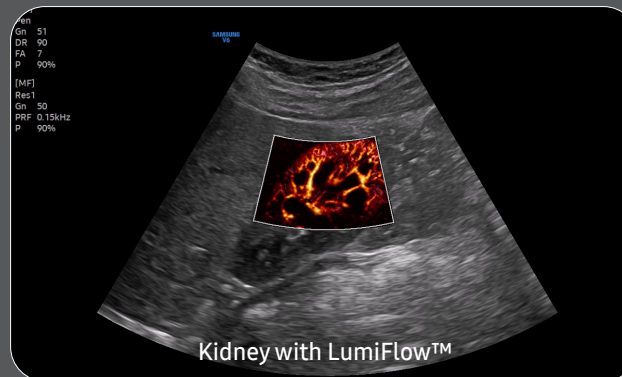

Kidney with LumiFlow™

Show blood flow in vessels in a 3D-like display

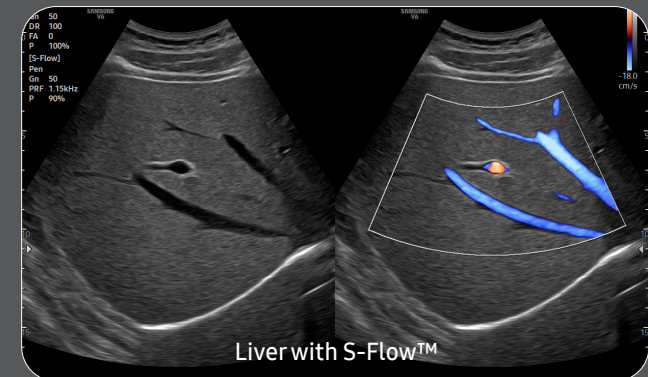

Liver with S-Flow™

Examine peripheral vessels with directional power Doppler

Reach new diagnostic confidence with comprehensive tools

Enhance your daily ultrasound diagnosis with the V6, a versatile solution created to efficiently support your clinical demands in general imaging. Benefit from our latest automation tools, which enable you to work with greater ease and achieve reliable results. Our aim is to assist you in prioritizing patient care, and the V6 stands as an ideal choice.




Display and quantify tissue stiffness in a non-invasive method

S-Shearwave Imaging™¹ allows the non-invasive assessment of stiff tissues in various applications. The color-coded elastogram, quantitative measurements, display options, and user-selectable ROI functions are useful for accurate diagnosis.



White paper

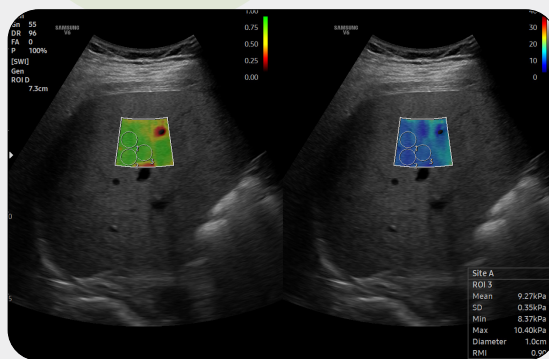


Quantitative measurement of liver fat with ultrasound signal

TAI™¹ (Tissue Attenuation Imaging) provides quantitative tissue attenuation measurement to assess steatotic liver changes.



White paper

TSI™¹ (Tissue Scatter distribution Imaging) provides quantitative tissue scatter distribution measurement to assess steatotic liver changes.

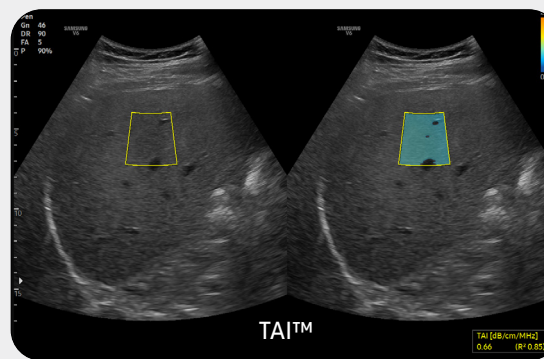

Hepato-renal index with automated ROI recommendation

HRI (Hepato Renal Index) is an index to quantify steatosis of a liver by comparing echogenicity between liver parenchyma and renal cortex. **EzHRI™¹** places 2 ROIs on the liver parenchyma and renal cortex and provides HRI ratio.



White paper



Analyze selected thyroid lesions and report thyroid assessment



S-Detect™^{1,3} for Thyroid analyzes selected lesions in the thyroid ultrasound study and shows the analysis data, provides standardized reporting based on the ATA, BTA, EU-TIRADS, and K-TIRADS* guidelines; and helps diagnosis with the streamlined workflow.

* ATA: American Thyroid Association

BTA: British Thyroid Association

EU-TIRADS: European Thyroid Imaging Reporting and Data System

K-TIRADS: Korean Thyroid Imaging Reporting and Data System

Quantify wall motion of the left ventricle

Strain+¹ is a quantitative tool for measuring global and segmental wall motion of the left ventricle (LV). Three standard LV views and a Bull's Eye are displayed in a quad screen for easy assessment of the LV function.



White paper

Detect and track nerves automatically with AI technology



NerveTrack™¹, a feature based on Deep Learning technology, detects and provides information of the location of the nerve area in real-time during ultrasound scanning.



White paper


Analyze selected breast lesions and report breast assessment



S-Detect™^{1,3} for Breast analyzes selected lesions in the breast ultrasound study and shows the analysis data, applies BI-RADS ATLAS* to provide standardized reporting; and helps diagnosis with the streamlined workflow.



White paper

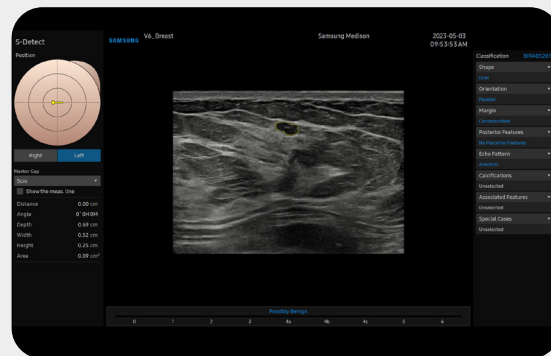

* Breast Imaging-Reporting and Data System, Atlas

It is a registered trademark of ACR and all rights reserved by ACR.


An automated reporting tool for heart diagnosis



HeartAssist™¹, a feature based on Deep Learning technology, provides automatic classification of ultrasound image into measurement views required for heart diagnosis and provides measurement results.



White paper

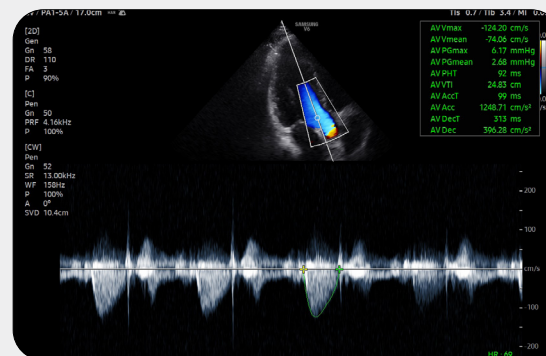

Easy calculation of the strain ratio between two ROIs

E-Strain™^{1,2} is designed to enable quick and easy calculation of the strain ratio between two regions of interest for day-to-day practice. Simply by setting the two targets, you can receive accurate, consistent results and make informed decisions in many types of diagnostic procedures.

Contrast Enhanced Ultrasound

CEUS+¹ is a contrast agent imaging technology. The micro-bubble contrast agent injected into the body through the vein or alike is subjected to perform nonlinear resonance due to stimulation of ultrasound energy.




White paper

Other features

AutoIMT+, AutoEF¹, StressEcho¹, ArterialAnalysis™¹, ElastoScan+™¹, NeedleMate+™¹, Panoramic+¹

Optimize workflow with precious time-saving tools

V6 is specifically designed to optimize the work efficiency of healthcare professionals. Notably through its remote accessibility, streamlined workflows, wider screen view for enhanced user experience, and its compact yet powerful design with battery capability, making it adaptable for diverse medical environments.

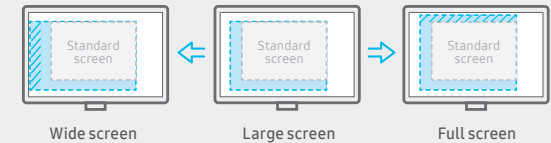


See images in expanded view

The ultrasound examination can be performed while viewing the images and cines that are expanded at various ratios according to the user preference.

Continue working even when AC power is temporarily unavailable

BatteryAssist™¹ provides battery power to the system, enabling users to perform scans when AC power is temporarily unavailable. It also allows the system to be moved without having to turn the power off and then back on.

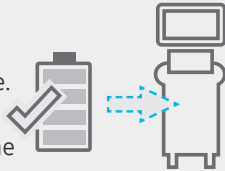

* The live scan time without power is about 3 times longer than the live scan time of the previous model.

Build predefined protocols to ensure every step is followed every time

EzExam+™¹ enables you to build or use a predefined protocol, and assign protocols for examinations that are regularly performed in the hospital in order to reduce the number of steps that you have to go through.




Customize frequently used functions on the touchscreen

TouchEdit, a customizable touchscreen, allows the user to move frequently used functions to the first page.



Compare previous and current exam in a side-by-side display

EzCompare™ automatically matches the image settings, annotations, and bodymarkers from the prior study.

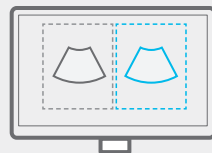


Select transducer and preset combinations in one click

QuickPreset allows the user to select the most common transducer and preset combinations in one click.




Assign functions to the buttons near the trackball

The buttons around the trackball can be customized for easy selection of commonly used functions.

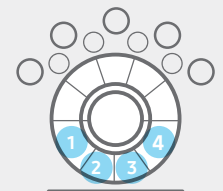

Save image data directly to USB memory

User can directly export image/cine with a USB device.

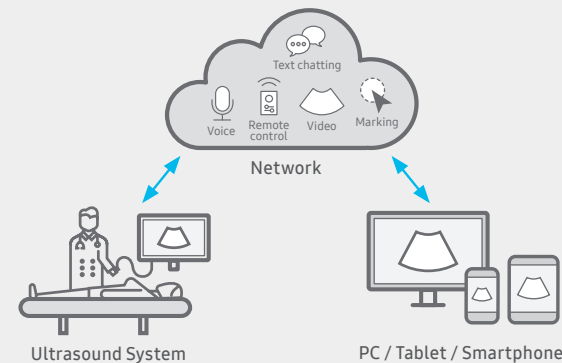

Real-time image sharing solution

SonoSync™ 1,5 is available in PC and smartphone, etc. as a real-time image share solution that allows communication for care guide and training between doctors and sonographers. In addition, voice chatting, text chatting and real-time marking functions are provided for better communication; and the MultiVue function is included that allows monitoring multiple ultrasound images on a single screen.



Learn more

Samsung healthcare cybersecurity

To address the emerging need for cybersecurity, Samsung provides a solution to support our customers by offering the tools to protect against cyberthreats that may compromise invaluable patient data and ultimately degrade the quality of care.



Learn more



Intrusion prevention



Access control




Data protection



Effective cooling system

An effective airflow system cools down the ultrasound system by constantly letting heat out and reducing fan noise.

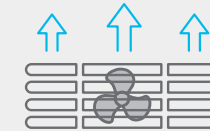

Recycled materials

Eco-friendly resin cover is applied to the air vent exterior cover, outlining Samsung's efforts towards a greener tomorrow.



Recycled materials

Eco Packaging

Reusable packaging composed of eco-friendly recycled paper. It is Samsung's commitment to achieving carbon-neutral of the earth and environment.



Learn more



Recycled materials



Eco Packaging

Comprehensive selection of transducers

Curved array transducers



Abdomen, Obstetrics, Gynecology, Pediatric, Musculoskeletal, Vascular, Urology



Abdomen, Obstetrics, Gynecology, Pediatric, Musculoskeletal, Vascular, Urology, Thoracic



Abdomen, Obstetrics, Gynecology, Pediatric, Musculoskeletal, Vascular, Urology, Thoracic



Abdomen, Pediatric, Vascular

Linear array transducers



Abdomen, Obstetrics, Gynecology, Pediatric, Musculoskeletal, Vascular, Urology, Thoracic



Abdomen, Pediatric, Musculoskeletal, Vascular, Small parts



Musculoskeletal, Intraoperative



Musculoskeletal, Pediatric, Vascular, Small parts

Phased array transducers



Cardiac, Vascular, Abdomen, Pediatric, TCD, Thoracic



Cardiac, Pediatric, Abdomen, Vascular, TCD



Cardiac, Pediatric, Abdomen, Vascular, TCD

Endocavity transducers



Obstetrics, Gynecology, Urology



Obstetrics, Gynecology, Urology



Obstetrics, Gynecology, Urology

Volume transducers



Abdomen, Obstetrics, Gynecology, Urology



Obstetrics, Gynecology, Urology

CW transducers



Cardiac, Vascular, TCD



Cardiac, Vascular, TCD

TEE transducers



Cardiac

* Ergonomic transducers

The new endocavity transducer supports natural grip by moving the max-width point to a more forward position and also increasing the length of the grip to allow balanced weight distribution.



Cleaning and disinfection guide

- * This product, features, options, and transducers may not be commercially available in some countries.
- * Sales and Shipments are effective only after the approval by the regulatory affairs. Please contact your local sales representative for further details.
- * This product is a medical device, please read the user manual carefully before use.
- * S-Vue Transducer™ is the name of Samsung's advanced transducer technology.

1. Optional feature which may require additional purchase.
2. Strain value for ElastoScan+™ is not applicable in the United States and Canada.
3. Recommendations about whether results are benign or malignant in S-Detect™ are not applicable in the United States.
4. SonoSync™ is an image sharing solution.

SAMSUNG MEDISON CO., LTD.

© 2023 Samsung Medison All Rights Reserved.
Samsung Medison reserves the right to modify the design, packaging, specifications, and features shown herein, without prior notice or obligation.

CE 0123